



To:

THE DIRECTOR GENERAL OF HEALTH SERVICES IN ALL GOVERNORATES
Commanding Officer, Armed Forces Hospital (Al Khoudh & Salalah)
Director General of Engineering Affairs, MOH
Director General of Royal Hospital
Director General of Khoula Hospital
Director General of Medical Supplies (MOH)
Director General of Pvt. Health Est. Affairs (to kindly arrange distribution to all Pvt. Hospitals)
Hospital Director (Al Nahda Hospital)
Hospital Director (Al Massara Hospital)
The Head of Medical Services in SQU Hospital
The Head of Medical Services in Royal Oman Police
The Head of Medical Services in Ministry of Defence
The Head of Medical Services in The Diwan
The Head of Medical Services in The Sultan's Special Force
The Head of Medical Services in Internal Security Services
The Head of Medical Services in Petroleum Development of Oman
The Head of Medical Services in LNG Oman
ALL PRIVATE PHARMACIES & DRUG STORES

After Compliments,

Please find attached our Circular No 195 dated 24/9/2023 Regarding NCMDR
FSN of Giraffe OmniBed and Giraffe OmniBed Carestation from (mfr: GE Healthcare).

Copy to:

- Director, Office of H.E. The Undersecretary for Health Affairs
- Director of Medical Device Control, DGPA&DC
- Director of Pharmacovigilance & Drug Information Dept, DGPA&DC
- Director of Drug Control Department, DGPA&DC
- Director of Pharmaceutical Licensing Department, DGPA&DC
- Director of Central Quality Control Lab., DGPA&DC
- Supdt. of Central Drug Information



Circular No. 195 / 2023

08 -03-1445 H

24 -09-2023

بتقنية بثقة
Moving Forward
with Confidence



Field Safety Notice of Giraffe OmniBed and Giraffe OmniBed Carestation from GE Healthcare.

| | |
|-----------------------|--|
| Source | NCMDR- National Center for Medical Devices Reporting- SFDA https://ncmdr.sfda.gov.sa/Secure/CA/CaViewRecall.aspx?caid=4&rid=19690 |
| Product | Giraffe OmniBed and Giraffe OmniBed Carestation. |
| Description | Incubator. |
| Manufacturer | GE Healthcare. |
| Local Agent | Muscat Pharmacy & Stores LLC. |
| The affected products | All Giraffe OmniBeds Giraffe OmniBed Carestation (2082844-001-XXX) [GTIN – 010084068211686221] Bedside panel FRU Kits: 5865384-01 to 5865384-34, 5865385-01 to 5865385-34, 5865386-01 to 5865386-34 and 5865387-01 to 5865387-34. |
| Reason | Potential an incorrect type of secondary latch could have been installed on the northeast (NE) and/or northwest (NW) bedside panels for certain devices. |
| Action | 1. Please inspect the device bedside panels as instructed in Appendix A of the attachment. If the device has the correct latches, continue to use the device. If the device does not have the correct latches, follow the instructions in Appendix A prior to clinical use. 2. Contact the local agent for remedial action. |
| Product Image | |
| comments | Healthcare professionals are encouraged to report any adverse events Suspected to be associated with the above device or any other medical device to Department of Medical Device Control through the E-mail: Med-device@moh.gov.om |

Dr. Mohammed Hamdan Al Rubaie

Director General



URGENT FIELD SAFETY NOTICE



Date of Letter Deployment

GE HealthCare Ref. # 32089


To: Director of Biomedical Engineering/Environmental Services
Director of Neonatology/Labor and Delivery/ Nurse Manager
Hospital Administrator

RE: **Giraffe OmniBed and Giraffe OmniBed Carestation potential incorrect secondary latch installed on certain bedside panels.**

Safety Issue

GE HealthCare has become aware of the potential that an incorrect type of secondary latch could have been installed on the northeast (NE) and/or northwest (NW) bedside panels for certain Giraffe OmniBed and Giraffe OmniBed Carestation devices (see figure 1).

Figure 1. Bedside Panels showing Northeast (NE) and Northwest (NW) latches

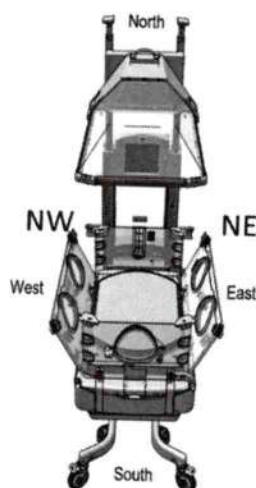

The Giraffe OmniBed and Giraffe OmniBed Carestation devices have primary and secondary latch mechanisms.

The primary latches are present and operate correctly.

The secondary latch mechanism serves two functions: (1) acts as a "catch" and reduces the potential for the bedside panel to fall if the user does not latch the primary latch mechanism correctly and (2) keeps the north wall in place when the canopy is raised.

For devices that have the incorrect type of secondary latch installed, the north wall can potentially become disengaged when the canopy is raised. However, the device

is designed such that the panels will stay closed and function as intended, even if only the primary latches on the south end are properly engaged.

There has been no customer complaints of this issue and no patient injuries.

**Actions
to be
taken by
Customer
/User**

Inspect the device bedside panels as instructed in **Appendix A**. If the device has the correct latches, continue to use the device. If the device does not have the correct latches, follow the instructions in Appendix A prior to clinical use.

Please complete and return the attached acknowledgement form to MIC.FMI32089@ge.com.

**Affected
Product
Details**

All Giraffe OmniBeds*

Giraffe OmniBed Carestation (2082844-001-XXX) [GTIN – 010084068211686221]
Bedside panel FRU Kits: 5865384-01 to 5865384-34, 5865385-01 to 5865385-34,
5865386-01 to 5865386-34 and 5865387-01 to 5865387-34

*NOTE: Some products were shipped prior to the implementation of UDI and may not contain a Global Trade Item Number (GTIN)

INTENDED USE:

The Giraffe OmniBed Carestation is a combination of an infant incubator and an infant warmer. The device can be operated as an incubator or as a warmer and can transition from one mode to the other on user's demand. It cannot be operated in both modes at the same time. Incubators and warmers provide heat in a controlled manner to neonates who are unable to thermo-regulate based on their own physiology. Incubators provide an enclosed, temperature-controlled environment and warmers provide infrared heat in an open environment. They may also be used for short periods of time to facilitate the neonate's transition from the uterus to the external environment. This device may incorporate a Servo Controlled Oxygen Delivery System. This is indicated to provide stable oxygen concentration within the infant compartment at the value set by the operator (21-65%).

**Product
Correction**

GE HealthCare will replace all affected latches at no cost to you.
A GE HealthCare representative will contact you to arrange for replacing any incorrect latches.
8004292222 SaudiArabiaServiceCenter@ge.com

**Contact
Information**

If you have any questions or concerns regarding this notification, please contact GE HealthCare Service or your local Service Representative.

Please be assured that maintaining a high level of safety and quality is our highest priority. If you have any questions, please contact us immediately per the contact information above.

Sincerely,



Laila Gurney
Chief Quality & Regulatory Officer
GE HealthCare




Scott Kelley
Chief Medical & Safety Officer
GE HealthCare

Appendix A

Northwest (NW) and Northeast (NE) bedside panel latch inspection

1. Visually inspect the NW and NE bedside panel latch assemblies on all Giraffe OmniBeds and Giraffe OmniBed Carestations to ensure the correct latches are present (See Figure 2).

Figure 2. North end bedside panel latch locations

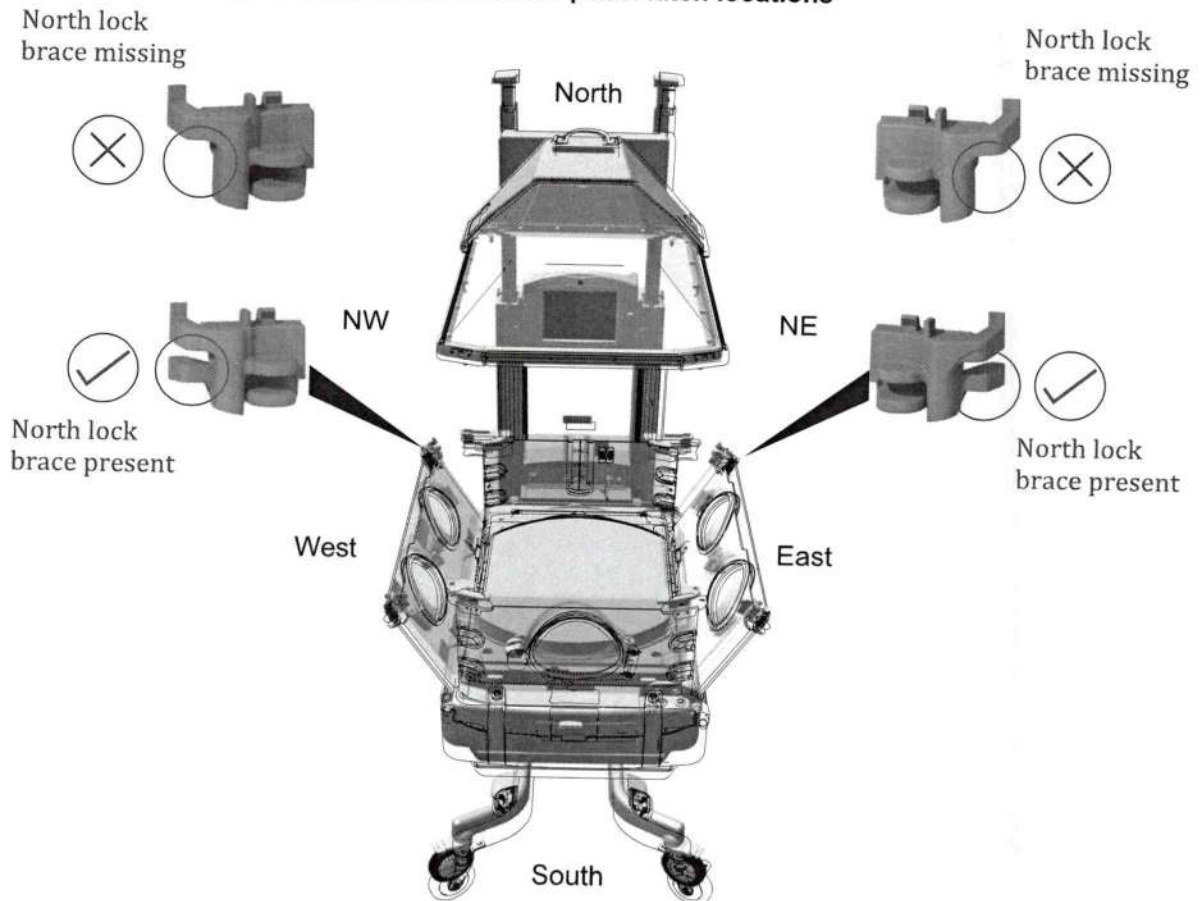

2. If an incorrect latch is found quarantine the device.
3. Inspect any field replaceable bedside panel stock at your facility per Figure 2 for the parts #'s listed below to ensure the correct latches are present and quarantine any field replaceable bedside panels with incorrect latches.
 - 5865384-01 to 5865384-34,
 - 5865385-01 to 5865385-34,
 - 5865386-01 to 5865386-34 and
 - 5865387-01 to 5865387-34
4. Record the results of the inspections on the provided Medical Device Notification Acknowledgement Response Form and send this form to GE HealthCare at MIC.FMI32089@ge.com
5. GE HealthCare will contact you to arrange for shipment of replacement latches. These will be provided free of charge to you.
6. Once correct latches are received, replace the incorrect latches with the correct latches and dispose of the incorrect latches.

If you need help with the above, please contact GE HealthCare Service at 1-800-437-1171.

Please provide the name of the individual with responsibility who completed this form.

Signature:

*Printed Name:

*Title:

*Date (DD/MM/YYYY):

*Indicates Mandatory Fields

Please return completed form by scanning or taking a photo of the completed form and email to: (MIC.FMI32089@ge.com)



**MEDICAL DEVICE NOTIFICATION ACKNOWLEDGEMENT
RESPONSE REQUIRED**

Please complete this form and return it to GE HealthCare promptly upon receipt and no later than 30 days from receipt. This will confirm receipt and understanding of the Medical Device Correction Notice. GE HealthCare needs this completed form to process shipment of replacement latches to you. These replacement latches will be provided free of charge.

*Customer/Consignee
Name: _____

Street Address: _____

City/State/ZIP/Country: _____

*Customer Email Address: _____

*Customer Phone Number: _____

We acknowledge receipt and understanding of the Medical Device Correction Notice and have executed the instructions as provided in this notification. We have inspected all of our **Giraffe OmniBed and Giraffe OmniBed Carestations and Field Replacement Stock** and provided the results below:

| Giraffe OmniBed and Giraffe OmniBed Carestation Serial Number or Field Replacement Stock Part Number | Latches are <u>correct</u> Yes or No |
|--|--|
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |
| | |

Please list every serial number and part number individually in the table above. If you need additional space to include serial numbers or part numbers, please copy this page and provide as many pages as necessary to include all devices.